



Ultrasonund Pocket DOPPLER S

PD Series

Model
PD-100Pro



Model
PD-1000Pro



Obstetrics and Gynecology

OB and Gynecology

Ultrasonic Pocket DOPPLERS

Superior Sensitiveness and High Fidelity Sound

Advanced® PD-100 Pro and PD-1000 Pro Ultrasonic Pocket Dopplers high performance dopplers with FHR digital LCD display with backlight and audio playback. All these models have audio output, which can be connected with earphone or recorder. As well, you can operate the unit with 4MHz or 8MHz probe for vascular diagnose.

Model PD-100Pro Ultrasonic Doppler



Key Features

- Easy to use
- Compact design
- Large backlight LCD display
- Mini USB Probe Socket
- Probe Detection
- Probe Identification
- Audio play & adjustable volume
- Earphone Socket
- Interchangeable waterproof probes
- Supports 2MHz, 3MHz, 4MHz, & 8MHz probes
- Automatic and manual counting of FHR (Fetal Heart Rate).
- Automatic shutdown
- Built-in Audio Recording & Playing
- Power by Alkaline Battery AA LR6 1.5 V.
- Power by NI-MH Battery Pack
- AC Battery Charger
- Battery Status Display
- Low Battery Indication
- Meets ISO 13485 Quality Standards

Model PD-1000Pro Ultrasonic Doppler



2 YEARS
warranty

Technical Specifications

Physical Specifications	Dimensions Weight	34mm Depth x 89mm Width x 141mm Height <300 g (including the battery)
Display	Resolution	45mm x 25mm LCD display
Probe	Nominal Frequency Specifications Cable length Weight Sensitivity Audio Output Power Protection Code	2.0 MHz Obstetrical Probe 3.0 MHz Obstetrical Probe 4.0 MHz Vascular Probe 8.0 MHz Vascular Probe 88mm (diameter) × 35mm (thickness) 2.5 meters 100 g 10 weeks gestation (3 MHz) 1 W IPX8 Water Ingress, indicating this probe does not get soaked under water within 1 meter deep for five hours
FHR Performance	Measuring Range Resolution Accuracy	50 bpm ~ 210 bpm 1 bpm ± 3 bpm
Recording and Playing	Audio Sampling Recording Length	Frequency 4KHz 240s
Control	White Backlight Auto Shut down	Two Brightness Adjustable: OFF, ON 1 minute after no signal or operation, auto shut down
Power Supply	Operating Voltage Frequency Alkaline Battery Rechargeable Battery Nominal Voltage Working Time Charge Time	100~ 240V 50/60 Hz (AA LR6 1.5 V) NI-MH battery 1800 mAh (optional) 2.4 VDC 8 hr 4 hr
Environment	Working Temperature Working Humidity Atmospheric Pressure	+5°C ~ +40°C (+41°F ~ +104°F) 25% ~ 80% (non-condensing) 860 hPa ~ 1060 hPa
Transport Storage	Temperature Humidity Atmospheric Pressure	-20°C ~ +55°C (-4°F ~ +131°F) 25% ~ 93% (non-condensing) 700 hPa ~ 1060 hPa
Safety	Anti-electric Shock Type Anti-electric Shock Degree	Internally powered equipment Type B equipment
PD-100 Pro	Standard Configuration	2MHz Waterproof Interchangeable Probe, AA battery, Carrying Case.
PD-1000 Pro	Standard Configuration	2MHz Waterproof Interchangeable Probe, Rechargeable battery, AC charger, Carrying Case.
Optional	Configuration	3MHz, 4MHz, 8MHz Probe, Earphones.

2018 Advanced Instrumentations Inc., is a U.S.A registered company – All rights reserved.

All functionality, features, specifications and other product information provided in this document including, but not limited to, the benefits, design, pricing, components, performance, availability, and capabilities of the product are subject to change without notice or obligation. Advanced Instrumentations reserves the right to make changes to this document and the product described herein, at any time, without obligation on Advanced Instrumentations to provide notification of such change. Actual description and specification of the product in this document may be different. Images shown here are for representational purpose only, actual may vary.


Advanced and Advanced Instrumentations trademarks and logos shown are property of Advanced Instrumentations Inc.

Success Through Quality/Since 1988

Advanced Instrumentations Inc.
Success Through Quality,
a Company You Can Trust

Advanced Instrumentations manufactures leading medical technology equipment in the areas of anesthesia, cardiology, operating room, gynecology and obstetrics, IV therapy, patient monitors, hospital furniture, neonatology and ultrasound. We deliver to the healthcare industry the highest-quality standards, reliability, and patient safety in all our products through effective, and rigorous testing procedures by our own department of Biomedical Engineering in the United States. All of our equipment comes with 2 years warranty and excellent post-sale support services.

Advanced Instrumentations Inc. Complies with the requirements and are audited by the US Food and Drug Administration (FDA) an entity of the health and Human Services of the United States of America. These certifications are the result of dedication and commitment to excellence in our products and services.

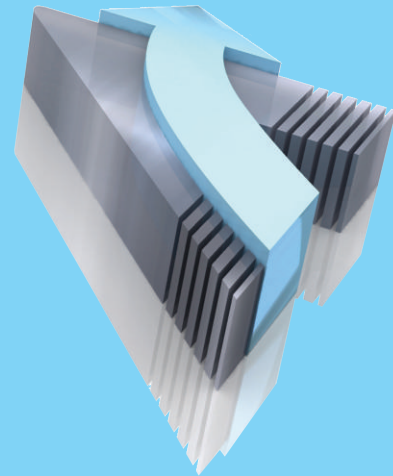

6800 N.W. 77 Court,
Miami, FL 33166
U.S.A.

Phone: 305-477-6331
Fax: 305-477-5351

For additional information visit us at:
www.advanced-inst.com